


Triton I

DENTAL TREATMENT UNITS

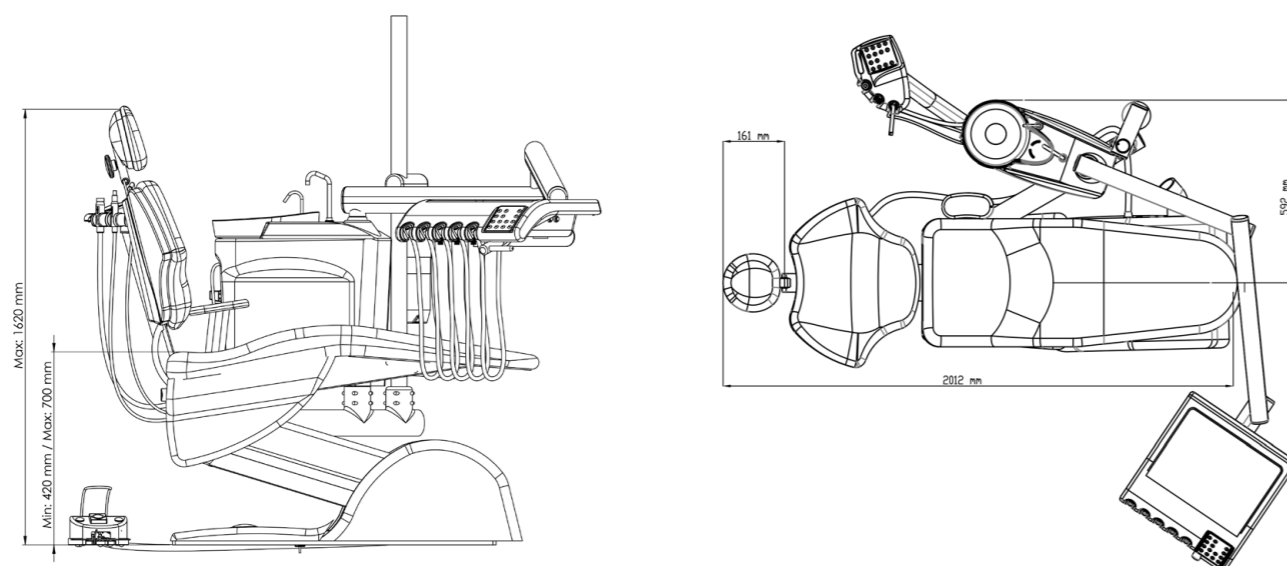

ELECTRICAL DATA

Input voltage	220/240 V AC
	Frequency 50/60 Hz
Energy consumption	100 to 900 VA
Electrical fuses	10 Amper
EU Certification	CE 0653


OPTIONAL EQUIPMENT

Right armrest	Available cables for digital external connections
Water supply bottle 2L	Assistant panel 3-way hydrosyrigne
Amalgam separator	

ENVIRONMENTAL CONDITIONS

Temperature	+10 to +40 °C
Relative humidity	30 to 75%
Air pressure	700 hPa to 1,060 hPa
Dental chair maximum load capacity	250 kg
Assistant tray maximum load capacity	5 kg


THE COLORS

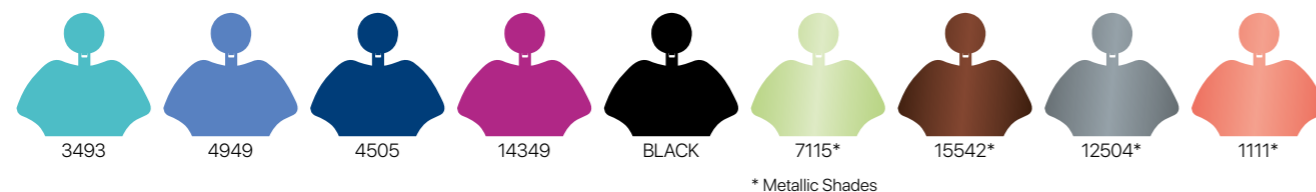

FACTORY
 10km. Old Highway
 Thessaloniki / Thessaloniki - Kilkis
 P.O. Box 142 GR 57008
 T: (+30) 2310 780 692 / (+30) 2310 780 792
 F: (+30) 2310 783 692
 E: info@neodent.eu

FOLLOW US ON



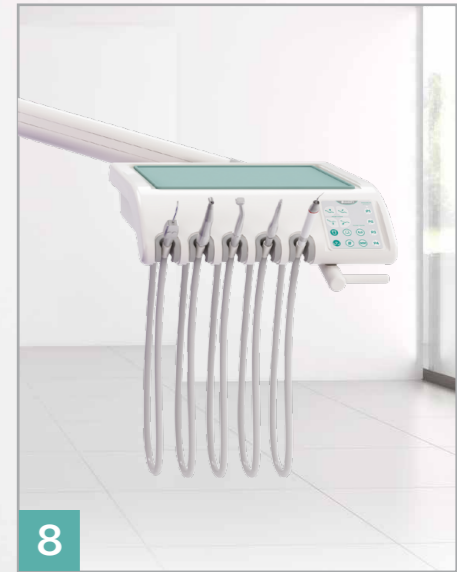

- 1** Adjustable Headrest, double-articulated, also suitable for children
- 2** Specially designed backrest; its clever shape – broader on top, narrower at the bottom- gives the feeling of robustness and comfort to the patient but at the same time allows the dentist to approach close to the chair.
- 3** Foot pedal, easy to use and move due to its metallic arch, offers multiple functions such as progressive instrument activation for the micromotor, chip-blower (optional), patient chair movement, last treatment position and hands-free operation for the dental light
- 4** Electromechanic system through the new, improved LINAK motors (Denmark), noiseless and smooth function. Safety Feature onto the chair base: in the presence of an obstacle the dental chair will automatically stop its movement
- 5** Synthetic upholstery, extremely durable, seamless, suitable for medical use, in a wide variety of colors
- 6** Nurse console with 2 aspiration hoses (diameter 11 and 17 mm) and 1 three-way syringe (optional), control panel for all functions as on the doctor's element.

- 7** Spittoon bowl, made of white porcelain, removable for cleaning and disinfecting. Its special curves ensure perfect sealing from water leakage.
- 8** Doctor's element in 3 different types: hanging cords, whip arm or mobile cart. Five -5- freely configurable outlets. Silicone cords, smooth and suitable for disinfection. Wide selection of rotating instruments: such as turbines, micromotors, piezo-electric scalers etc.
- 9** Soft touch control panel with buttons for each function: movement of the seat, 3 programmable memories, P1 rinsing position, water to bowl, water to cup, on/off water to instruments, auxiliary button.
- 10** Extremely flexible cords at the whip arm edition that can be secured through a locking system into the desired working position for optimum comfort to the operator's hands.
- 11** Instrument tray of dimensions 36 X 28 cm as additional storage space. The autoclavable rubber mat protects the working surfaces from undesired scratches.